

SC760-US Suction Controller



The innovative lightweight design and color coded body provide excellent performance with uncompromising safety

- ◆ Lightweight construction
- ◆ Color coded models
- ◆ Superior, integrated pipeline overflow protection with viral filter option
- ◆ Separate on/off switch
- ◆ Precise linear action vacuum adjustment and display
- ◆ MRI Compatible*

Excellence and Innovation


Additional Features

- ◆ Venturi injector suction unit also available (oxygen or air)
- ◆ Intermittent suction unit also available
- ◆ 12 month warranty
- ◆ Coordinates with our oxygen therapy range
- ◆ Pipeline overflow protection
 - BactiTrap® bacterial filter and mechanical float overflow protection (fitted as standard)
 - HydroTrap, bacterial/viral filter and hydrophobic float overflow protection (order separately and fit as a direct replacement for BactiTrap®)


Ordering Information

SC760-US Suction Controller

<i>Part No.</i>	<i>Description</i>
310128	High Suction with 1/8" NPT Universal thread, and BactiTrap®
310129	Low Suction with 1/8" NPT Universal thread, and BactiTrap®
4830-014	BactiTrap®, bacterial filter and over-flow protection
9970-001	BactiTrap®, box of 48
310025	HydroTrap, bacterial/viral filter and over-flow protection
310027	HydroTrap, box of 48

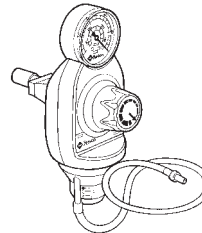


310128

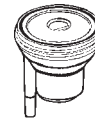

4830-014

AV500 Injector Suction Unit

<i>Part No.</i>	<i>Description</i>
312040	AV500 Injector suction unit with 1/4" inlet connection and MAK 500 receiver jar

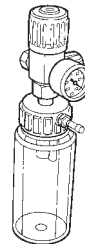

312040

*MRI Compatibility

MRI compatibility does not infer zero magnetism but an acceptable level of magnetism within a specified magnetic field.

The High Vacuum Suction Controller, when used within 1000 Gauss of magnetic exposure, functions within specification, does not represent a projectile danger and does not distort the MRI image.

Distributed by:

Suction Control Technology



Penlon Limited, Abingdon Science Park, Barton Lane, Abingdon, OX14 3PH UK

Sales Office

Tel: +44 1235 547001
 Fax: +44 1235 547021
 E-mail: export@penlon.co.uk

Penlon Limited reserves the right to change technical specification without prior notice.

© 2004 Penlon Limited



Penlon Limited is a Member of the InterMed Group